

# ***Aculus-AM-320***

Digital Honeycomb Laser Hood



## **User's Guide**

English

Getting started	1
Using Aculas-AM-320	4
Tips and troubleshooting	5
Learning more, service and support	6
Safety and cleaning	7


# Getting started

Read on to find out how to get started low level laser therapy (LLLT) with Aculas-AM-320

---

## 1 Connect Aculas-AM-320

Connect your Aculas-AM-320 power cable to a working electric outlet.

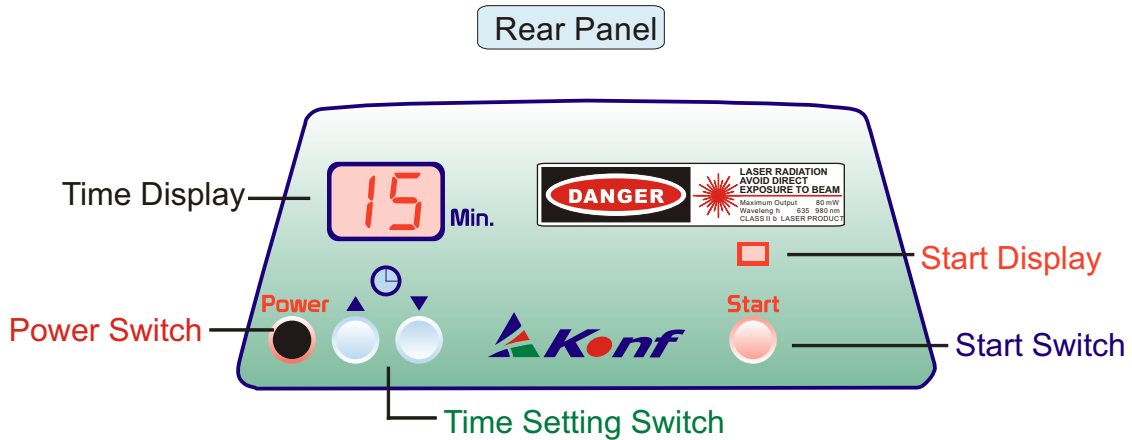

## 2 Power ON display








Turn on Aculas-AM-320's power switch and the power ON display will indicate the status.

Front Panel

Power ON Display



### 3 Operation steps and LED display

Operation			LED Display
STEP	ACTION	PUSH	
1	Power ON		
2	Start		
3	Time Setting	 	
Auto Time Setting: 15 Minutes Manual Time Setting: 1 ~ 60 Minutes			

# Using Aculas-AM-320

Aculas-AM-320 can be used for treatment of diseases and conditions, follow the recommended dosage and treatment time to get the best efficacy. You can also find new information about using Aculas-AM-320 on the web [www.konftec.com](http://www.konftec.com).

---

## Feet Assembly



## Operate Extension Rod



## Warning



# Tips and troubleshooting

Read on for advice about using and troubleshooting Aculas-AM-320.

---

## **My Aculas-AM-320 won't turn on**

- ☑ Make sure Aculas-AM-320's power switch is on.
- ☑ If that doesn't work, check the connection from Aculas-AM-320 to a working electrical outlet.
- ☑ If that doesn't work, contact your local sales office.

## **My Aculas-AM-320 is continuously playing music**

- ☑ **Turn off the power switch** then reconnect the power cord.
- ☑ Repeat the operation steps:
  1. Turn on Aculas-AM-320's power switch.
  2. Press "Start" button to begin the therapy.
- ☑ If that doesn't work, contact your local sales office.

## **Learning more, service, and support**

There is a wealth of information about using Aculas-AM-320 on the web.

For the latest information on Aculas-AM-320 service and support information go to

[www.konftec.com](http://www.konftec.com).

---

### **Obtaining warranty service**

If the product appears to be damaged or does not function properly, please follow the advice in this booklet, and the online resources. If the unit still does not function, go to

[www.konftec.com](http://www.konftec.com) for instructions on how to obtain warranty service.

# Safety and cleaning

Read on to learn about using Aculas-AM-320 safely and cleaning Aculas-AM-320.

---

## Important safety instructions

When setting up and using your Aculas-AM-320, remember the following:

- ☑ Read user's guide carefully before you plug your Aculas-AM-320 power cord into a power outlet.
- ☑ Keep user's guide handy for reference by you and others.  
Follow all instructions and warnings dealing with your Aculas-AM-320.
- ☑ **Warning:** Electrical equipment may be hazardous if misused.  
Operation of this device, must always be supervised by an adult. Do not allow children access to the interior of any electrical product and do not permit them to handle any cables.

## Warning: Laser beam

- ☑ Avoid direct eye exposure to the laser beam.
- ☑ Because some Aculas-AM-320 lasers emit invisible infrared light, the usage of Aculas-AM-320 must proceed with appropriate care and caution. Invisible infrared laser light penetrating the eye may cause serious harm after extended exposure. Therefore, the use of the optional laser safety goggles for the patient and clinician is advised.



### **Using the power cord**

- ☑ The only way to shut off power completely is to disconnect it from the power source.
- ☑ When connecting or disconnecting your power adapter, always hold the power adapter by its sides. Keep fingers away from the metal part of the plug.

### **About operating and storage temperatures**

- ☑ Operate your Aculas-AM-320 in a place where the temperature is always between 0 and 40° C (32 to 105° F).
- ☑ Store your Aculas-AM-320 in a place where the temperature is always between -20 and 85° C (-4 to 185° F).

When you are using your Aculas-AM-320, it is normal for the case to get warm.

## **Avoid wet locations**

**Warning** To reduce the chance of shock or injury, do not use your Aculas-AM-320 in or near water or wet locations.

- ☑ Keep your Aculas-AM-320 and power adapter away from sources of liquids, such as drinks, wash basins, bathtubs, shower stalls, and so on.
- ☑ Protect your Aculas-AM-320 and the power adapter from direct sunlight and rain or other moisture.
- ☑ Take care not to spill any food or liquid on Aculas-AM-320 or its power adapter. If you do, unplug Aculas-AM-320 before cleaning up the spill.
- ☑ Depending on what you spilled and how much of it got into your equipment, you may have to send your equipment to Konftec for service. See “Learning more, service, and support.”

## **Do not make repairs yourself**

**Warning** Do not attempt to open your Aculas-AM-320, disassemble it. You run the risk of electric shock and voiding the limited warranty. No user-serviceable parts are inside. For service, see “Learning more, service, and support.”

## **Cleaning**

Follow these general rules when cleaning the outside of your Aculas-AM-320:

- ☑ Make sure your Aculas-AM-320 is unplugged.
- ☑ Use a damp, soft, lint-free cloth. Avoid getting moisture in openings.  
Don't use aerosol sprays, solvents, alcohol, or abrasives.

## **About handling**

Your Aculas-AM-320 may be damaged by improper storage or handling.

Be careful not to drop your Aculas-AM-320 when operating or transporting the device.



[www.konftec.com](http://www.konftec.com)

© 2008 Konftec Corporation. All rights reserved.  
Konf, the Konf logo, Aculas and Aculas-AM  
are trademarks of Konftec Corporation,  
registered in Taiwan and  
other countries.

UAM320\_81001E